


LiNA Bipolar Loop™

*Laparoscopic Supracervical Hysterectomy (LSH)
and subserosal pedunculated myomectomy*

Amputation in a few seconds



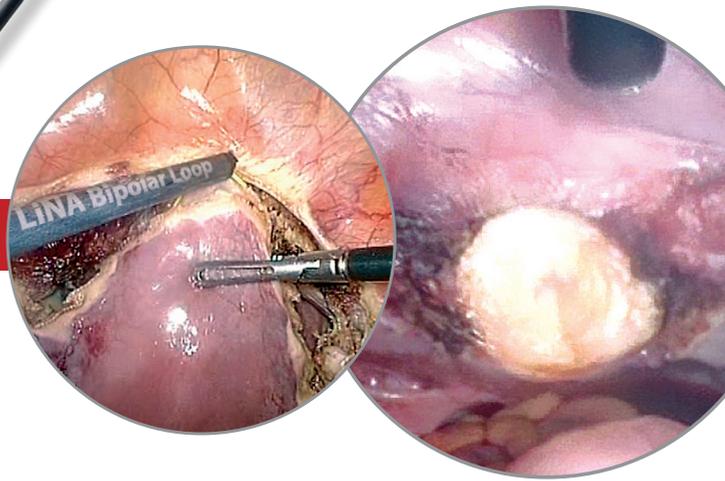

LiNA Bipolar Loop™

*Laparoscopic Supracervical Hysterectomy (LSH)
and subserosal pedunculated myomectomy*

Safety. Low energy, high effect.

Bipolar technology offers immediate cutting action and improved safety with only 55-70 watts. No grounding pad needed.

Visualization is aided by a protective ceramic tip, allowing for optimum accuracy in loop positioning.

The protective ceramic tip in conjunction with the yellow insulation isolates the cutting action to the tissue between the two electrodes for improved safety.

Fast. Clean cut in seconds.



Reduce your cervix amputation time by at least several minutes.

Cut and coagulation in one step leaves a clean hemostatic surface, which reduces the need for further coagulation of the cervical stump or peduncle.

Simple. No special generator required.

LiNA Bipolar Loop™ is compatible with standard electrosurgical units to assure quick and easy setup.

3 meter cable with standard 2-pin bipolar plug.

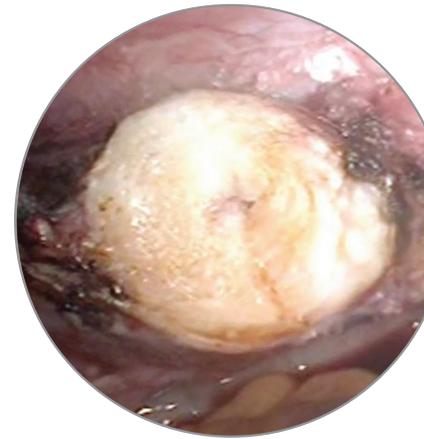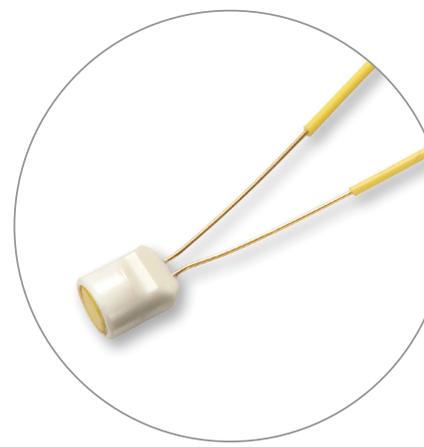

Ordering Information:

Product code	Description	UOM
BL-160	LiNA Bipolar Loop™ – Endoscopic Bipolar Loop for LSH & Myomectomy. Small Loop - 160mm x 80mm. 2-pin plug. Sterile, single use	Box of 6 ea
BL-200	LiNA Bipolar Loop™ – Endoscopic Bipolar Loop for LSH & Myomectomy. Large Loop - 200mm x 100mm. 2-pin plug. Sterile, single use	Box of 6 ea
BL-240	LiNA Bipolar Loop™ – Endoscopic Bipolar Loop for LSH & Myomectomy. X-Large Loop - 240mm x 120mm. 2-pin plug. Sterile, single use	Box of 6 ea

Distributor:

LiNA Medical USA

Norcross, GA 30093

Tel: +1 (855) 546-2633 | Fax: +1 (678) 380-0500 | info@linamed.com | www.linamed.com

©2025 LiNA Medical AG. All rights reserved. All products may not be available for sale in all geographies. Please contact your local LiNA Medical representative for information about products available in your area. This information is intended for healthcare professionals. If you are a patient seeking more information, please consult your healthcare provider. Before using any LiNA product, always refer to the IFU included in our packaging for complete indications for use, contraindications, warnings, precautions, and instructions.

LiNA is a registered trademark of LiNA Medical in EU, US and CN.



LiNA Medical AG, Platz 3, 6039 Root D4 | Switzerland
Phone: +41 78 617 56 56, info@lina-medical.com | www.linamed.com



LiNA Medical Polska Sp. z o.o., Rolna 8A, Sady, 62-080 Tarnowo Podgórze | Poland
Phone: +48 61 222 21 21, info.poland@lina-medical.com | www.lina-medical.pl

EN ISO 13485:2016
Quality Management System

