


# LiNA **Gold Loop**®

Laparoscopic Supracervical Hysterectomy (LSH)



Amputation in a few seconds

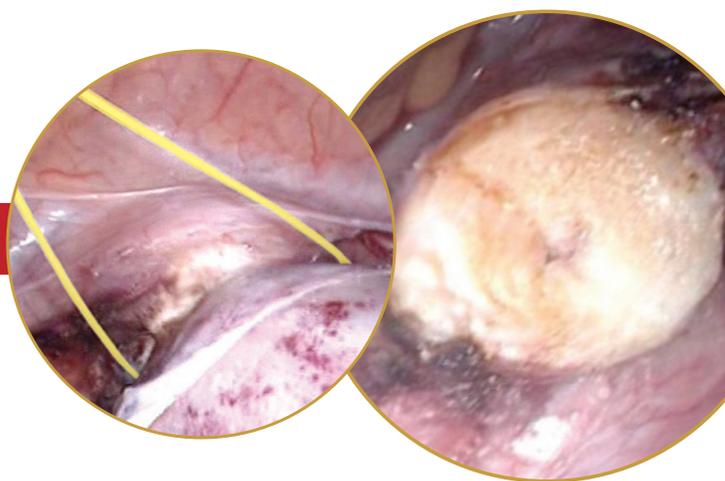

Innovation in Gynecology


# LiNA Gold Loop®

## Laparoscopic Supracervical Hysterectomy (LSH)

### Efficiency. Minimize procedure time by at least several minutes.

The LiNA Gold Loop® is a safe and effective monopolar loop for sectioning of the uterus during laparoscopic supracervical hysterectomy (LSH).


Cut and coagulate in one single step. Minimize the necessity for post coagulation of the surface by utilizing a blend mode (the recommended setting for the generator is: 100-130 W, Blend 1).

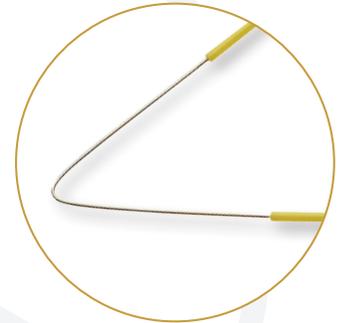

### Versatility. 3 sizes available and no special generator required.

The LiNA Gold Loop® is compatible with standard electrosurgical units to assure quick and easy setup.

3 meter cable with standard 8 mm male connector included.

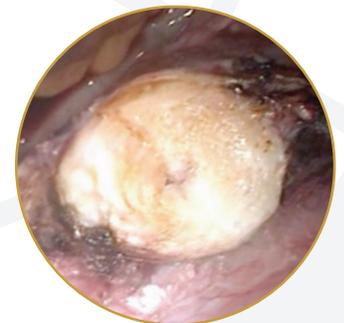

### Safety. Precise sectioning of the uterus.

Yellow insulation on the loop improves visibility on the monitor.

Provides a clean and coagulated surface.



**Important:**  
Please refer to the Instruction for Use for complete instructions, contraindications, warnings and precautions.

#### Ordering Information:

Product code	Description	UOM
EL-160-8	<b>LiNA Gold Loop®</b> – Endoscopic Monopolar Loop for LSH. Small Loop - 160 mm x 80 mm, 8 mm connector. Sterile, single use.	Box of 6 ea
EL-200-8	<b>LiNA Gold Loop®</b> – Endoscopic Monopolar Loop for LSH. Large Loop - 200 mm x 100 mm, 8 mm connector. Sterile, single use.	Box of 6 ea
EL-240-8	<b>LiNA Gold Loop®</b> – Endoscopic Monopolar Loop for LSH. X-Large Loop - 240 mm x 120 mm, 8 mm connector. Sterile, single use.	Box of 6 ea

**Distributor:**  
**LiNA Medical USA**  
Norcross, GA 30093  
Tel: +1 (855) 546-2633 | Fax: +1 (678) 380-0500 | [info@linamed.com](mailto:info@linamed.com) | [www.linamed.com](http://www.linamed.com)

©2025 LiNA Medical AG. All rights reserved. All products may not be available for sale in all geographies. Please contact your local LiNA Medical representative for information about products available in your area. This information is intended for healthcare professionals. If you are a patient seeking more information, please consult your healthcare provider. Before using any LiNA product, always refer to the IFU included in our packaging for complete indications for use, contraindications, warnings, precautions, and instructions.

LiNA is a registered trademark of LiNA Medical in EU, US and CN.

 LiNA Medical AG, Platz 3, 6039 Root D4 | Switzerland  
Phone: +41 78 617 56 56, [info@lina-medical.com](mailto:info@lina-medical.com) | [www.linamed.com](http://www.linamed.com)

 LiNA Medical Polska Sp. z o.o., Rolna 8A, Sady, 62-080 Tarnowo Podgórze | Poland  
Phone: +48 61 222 21 21, [info.poland@lina-medical.com](mailto:info.poland@lina-medical.com) | [www.lina-medical.pl](http://www.lina-medical.pl)

EN ISO 13485:2016  
Quality Management System

